

MediCap™ USB200

Medical Digital Video Recorder




The easy way to record video and images

The MediCap USB200 is the medical industry's easiest full-featured video recording solution. Your videos and still images are recorded directly to a USB flash drive or optional hard drive. The recorded videos are DVD-quality, but without all the hassles of DVD disks. The USB200 is so simple to use, you'll be recording videos and images just minutes after you unpack it. Now it's easy to step up to solid-state digital recording – and at a surprisingly affordable price.

Powerful and affordable...

- ▶ Works with endoscopes, arthroscopes, ultrasounds, or any video source with standard outputs
- ▶ Record an hour of video on the included flash drive or hundreds of hours on an optional external hard drive
- ▶ DVD-quality video, plus images in JPG, DICOM, etc.
- ▶ Collect patient information using any USB keyboard
- ▶ Record via button, footswitch, or trigger
- ▶ Network connection option
- ▶ Compact, rugged design fits easily on any cart
- ▶ No training required – a snap to connect and use

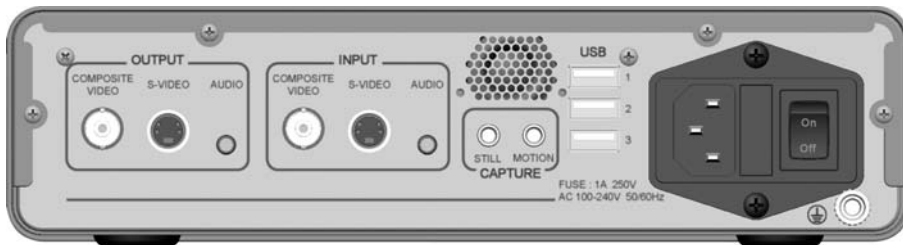

MediCap™ USB200 Specifications

Device Type	Medical Digital Video Recorder Records digital video and still images from analog video sources to USB flash drive or optional USB hard drive
Case	Rugged metal case with high-impact plastic faceplate
LCD Monitor	2 1/2", full color (for navigation only, not for diagnosis)
Speaker	200Hz-15kHz, 1W, (for audio confirmation only, not for diagnosis)
Buttons	Sealed membrane, fluid resistant
Video In/Out	PAL or NTSC (menu selectable) S-video or composite BNC connections
Audio In/Out	Standard 3.5mm stereo phone jacks
USB Support	USB2.0, high speed, 1 front connector, 3 rear connectors (all powered)
Networking	Ethernet T10/100 connection via Network Upgrade Kit (optional)
Footswitch	Hands-free capture Standard 3.5mm mini jack
Mechanical	Size: 240 x 200 x 63mm / 9.5 x 7.9 x 2.5" Weight: 1.7kg / 3lb 12oz (appx.)
Environmental	Storage: -40° to +85° C Operating: -20° to +40° C
Power	Power 110/220V 3-prong AC jack
Media	USB flash drives or external USB hard drive
Image Formats	JPEG, TIFF, PNG, DICOM
Image Size	1024 x 768, 800 x 600, 640 x 480 pixels
Video Format	MPEG2 (DVD-quality)
Audio Format	MPEG1 L2, 2 channels
Certifications	UL Medical, CE, FCC, FDA Class 1, RoHS
Kit Contents	MediCap USB200 USB flash drive Video cables Power cable User guide

Also Available...



MediCap USB170

All the power and features of the USB200, but for still images only, and with an even more affordable price.